Accelerate your teeth whitening with the use of light therapy. Our whitening device can be used to whiten your teeth faster and without the pain that commonly occurs with bleaching.

YOUR KIT INCLUDES:

1

BEFORE STARTING NOTE:

- For better results, we suggest having your teeth cleaned by a dental hygienist prior to using a whitening kit.
- If teeth become overly sensitive, stop treatment.
- The gel removes years of stains caused by coffee, tea, wine, smoking, etc.
- Avoid directly applying teeth whitening gel to gum area.
- Wipe off excess.
- If gel gets on your face or clothes, you should wash them without any delay.

2

CHARGING DEVICE SPECIFICATIONS:

Charging time - **2** hours
One charge can last **5-6** uses


STORAGE:

- Keep LED Light away from sunlight and keep in cool and dry area.
- The gel pens should be kept in the refrigerator.

instructions

- 1 Brush and floss teeth then rinse mouth prior to using mouthpiece.
- 2 While looking in the mirror, compare your teeth to shade guide paper to record which stage your teeth are in.
- 3 Rinse mouth piece under warm water and apply whitening gel by painting the gel to your upper and lower teeth avoiding lips, gums or cheeks.
 - * Do not apply too much gel, more gel will not get you faster or better results.
 - * If gel comes into contact with lips, gums or cheeks, you may feel a little burning sensation, rinse with water and the feeling will go away within 24 hours.
- 4 Place mouthpiece in mouth, bite down, take deep breaths. Breathe through your nose.
- 5 Turn on the light by tapping on the power button. The whitening device will start on STAGE 1 of the whitening process with the blue light for 15 minutes.
 - * The blue light is used to activate the gel quickly without any heat or ultraviolet light, working to weaken yellow stains.
- 6 After stage one is complete, mouthpiece will turn off. Tap power button again and STAGE 2 will begin with the blue and red light for 10 minutes.
 - * Red light is used to stimulate cellular regeneration which helps heal minor mouth irritations.
- 7 Once both stages are complete, the whitening device will turn off by itself, remove from mouth and rinse mouthpiece under warm water.
- 8 Compare teeth to shade guide paper again to see results of how many stages your teeth improved.
 - * Do not eat or drink 30 minutes after you've finished treatment.

START TO WHITEN

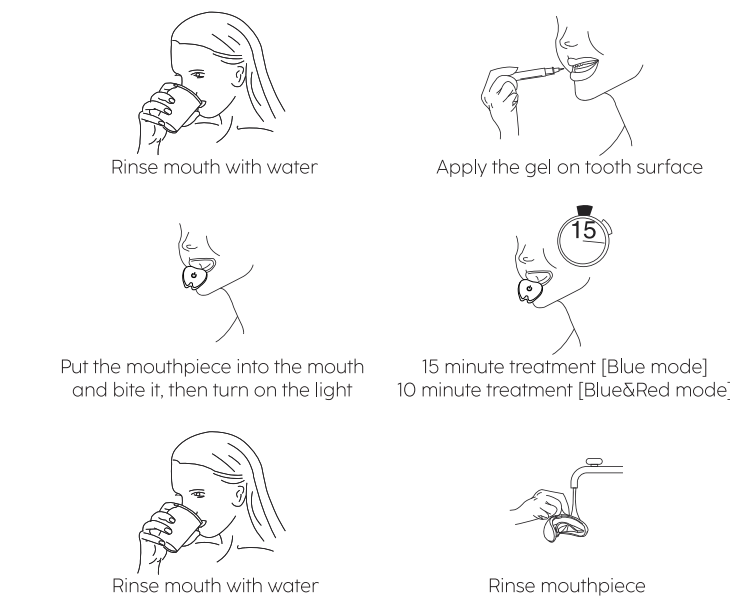Gel Pen Ingredients:
Carbamide Peroxide, Glycerol, Propylene Glycol, Deionized Water, Carbomer, Carboxymethyl, Poly Vinylpyrrolidone, Sodium Hydroxide, Menthol.

4

TERMS AND CONDITIONS

Do not use if you are under the age of 18. Do not use if you have overly sensitive teeth. The product should only be used to whiten your teeth. Do not use if your teeth are decaying or loose. Whitening product will not whiten any existing dental restoration such as fillings, crowns, veneers, bridges. Do not use if you are pregnant or lactating. Do not use if allergic to the ingredients listed above. Keep out of reach of children.

5

Bringing you white smiles through the power of light technology

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